




No More Flaking

CHEMVAC™ TEFLON® STIRRER BEARING

CG-2077-G STIRRER BEARING, CHEMVAC™, TEFLON®
NEW Chemvac™ Teflon stirrer bearing provides a VACUUM TIGHT FLAKE-FREE SEAL. The bearing consists of a Teflon standard taper body, Glass filled polypropylene screw cap, Peek Compression spring, Teflon/Peek sleeve and a glass filled polypropylene loosening nut. The polymer blend sleeve provides a long lasting seal for use up to 500rpm. Loosening nut above the standard taper inner joint allows easy and quick removal. *NOTE: Teflon stirrer shafts are not recommended for use with this bearing.*

Features:

- FLAKE-FREE
- Chemically Resistant
- Vacuum (~5 Torr) And Pressure (~14 PSI)
- Polymer Blend Sleeve Provides Excellent Wear Resistance
- For Use With Polished Glass Or Stainless Steel Stirrer Shafts

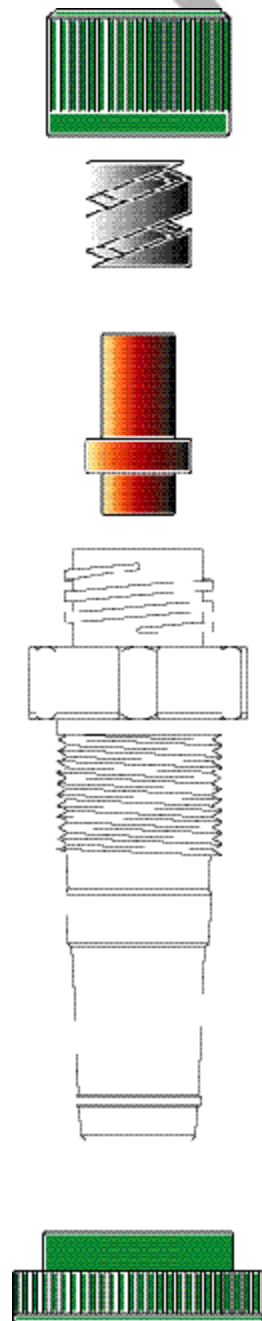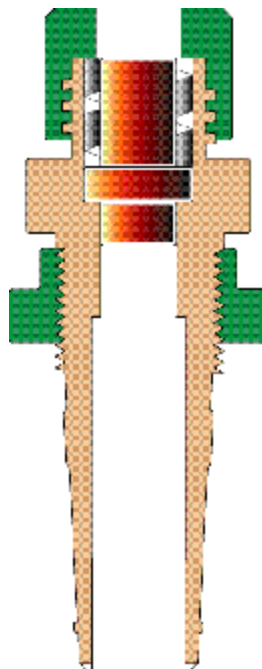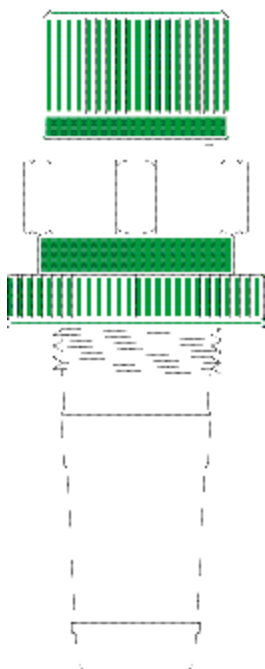

<u>Part Number</u>	<u>Description</u>	<u>Price</u>
CG-2077-G-01	24/40 Teflon bearing, 10mm	112.50*
CG-2077-G-02	29/42 Teflon Bearing, 10mm	118.20*
CG-2077-G-04	45/50 Teflon Bearing, 10mm	133.00*

Replacement Components:

<u>Part Number</u>	<u>Description</u>	<u>Price</u>
CG-2077-G-10	Replacement Bearing Sleeve, 10mm	18.00*
CG-2077-G-11	Replacement Peek Compression Spring, 10mm	19.70*

*NOT SUBJECT TO DISCOUNT

CHEMGLASS
Scientific Apparatus Since 1946

3861 North Mill Road • Vineland NJ 08360
U.S. and Canada phone: 1-800-843-1794
Fax: 1-800-922-4361

customer-service@chemglass.com • www.chemglass.com