

IMPROVED SUPPORT CLAMP


CG-1947-A SUPPORT CLAMP, STAINLESS STEEL, IMPROVED

Similar to CG-1947, support clamp is designed to fit just below the flange of all Chemglass Reactors and permits Reactor lid to be detached without removing the Reactor body from the support clamp. Improved clamp is hinged on one side to swing open and allow complete access to install/remove reactor while clamp is still attached to lattice or benchtop support stand. Removable tightening knob holds Reactor secure. Clamps are constructed of stainless steel with a 1/2" OD X 9.5" long extension rod for connection to lattice and benchtop support stands. The inner portion of the clamp is lined with a high temperature red silicone rubber to prevent direct contact between the stainless steel and glass Reactor.

PLEASE NOTE: IN ADDITION TO USING A CG-1947-A CLAMP, ALL REACTORS SHOULD BE SUPPORTED FROM THE BOTTOM.

Part Number	Fits Reactor with Flange Size (mm)	Price
CG-1947-A60	60	98.00*
CG-1947-A80	80	105.00*
CG-1947-A100	100	110.00*
CG-1947-A150	150	115.00*

*not subject to discount

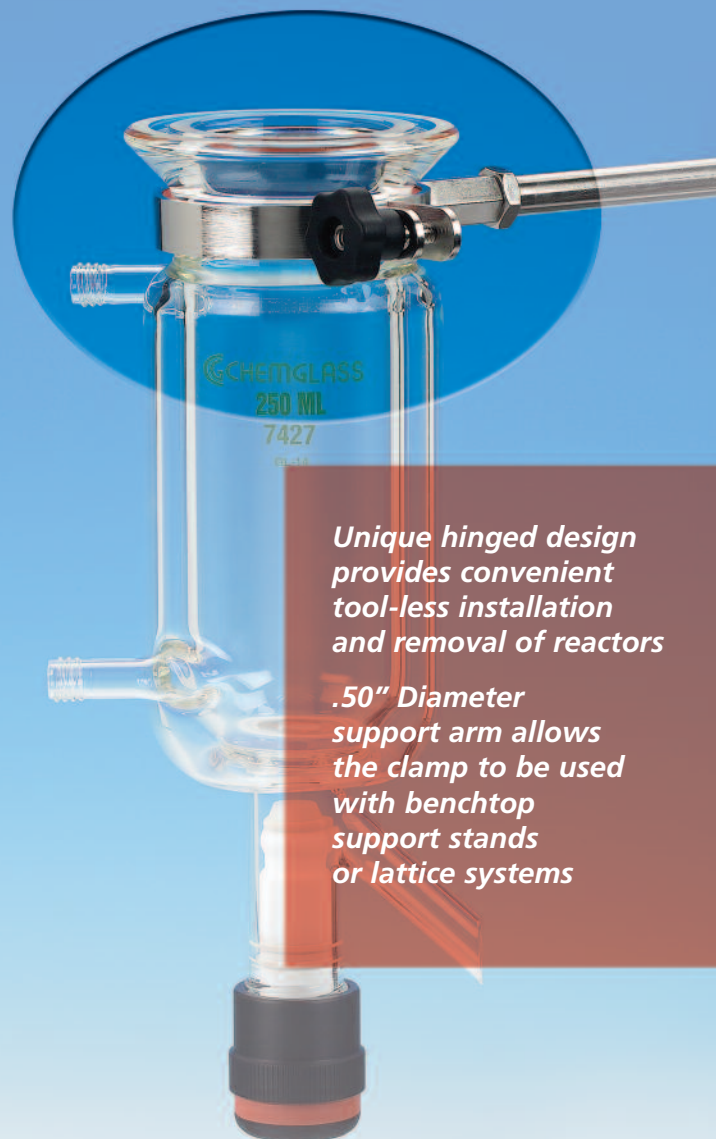

Unique hinged design provides convenient tool-less installation and removal of reactors

.50" Diameter support arm allows the clamp to be used with benchtop support stands or lattice systems



CG-1947-A